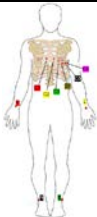

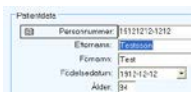
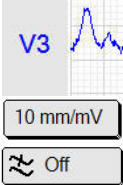




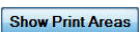




Quick start guide Long ECG *EC Sense*

This quick start guide gives a brief explanation on how to use *EC Sense* for acquiring a Long ECG.

Icon	Fn	Instruction
		1. Prepare the patient and connect all electrodes according to the desired lead electrode placement. Check that the patient cable is connected to the acquisition unit.
	F2	2. Click on the icon for entering patient data.
		3. Enter the patient information in the dialogue box. If the PC is connected to a central storage system the information may be filled in automatically.
		4. Select desired leads, gain and filter settings. Check the quality of the ECG signals.
	F9	5. Click on the start icon to start the ECG recording.
	F11	6. Enter blood pressure and other values, if desired.
	F10	7. Mark the events that should be included in the report.
	F9	8. Stop the recording by clicking on the red square icon. Disconnect the patient.
		9. In the report mode you can decide if you wish to add events and include them for printout before the report is printed and/or saved.
	F4	10. Click on the printer icon to printout the entire recording (or use <i>Print page</i> to print specific pages).
	F5	11. Save the recording on the local PC or in the central storage system.
Read the Operator's manual for complete instructions!		



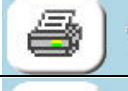
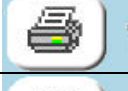

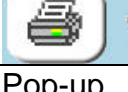






CARDIOLEX

Vretenvägen 6
171 54 Solna
Sweden

Version 4.0

Commands in Long ECG mode

The summary below shows the available functions in the *Long ECG* mode.

Function	Icon	Menu selection	Fn key	Comment
Patient Entry		File → Patient Info	F2	
Stop examination		File → Exit		No ECG data will be saved
Real Time ECG on printer		File → Man Print	F5	
Rhythm ECG on printer		File → Rhythm Print	F6	
Print a Status Report		File → Status Report	F4	
Print Latest Event/Arrhythmia			Alt+F4	Event/Arrhythmia must first be added
Print a ECG interval	Pop-up in rhythm	Print from/to...		Print from time lock to clicking point
Start recording		File → Long ECG	F9	Alter between Start/Stop
Stop recording			F9	ECG is saved and reports are opened
Mark an event on ECG		Actions → Mark event	F10	
Enter trend data		Actions → Enter measurements	F11	
Take the blood pressure and enter the BP values		Actions → Enter BP	F12	
Select or deselect tremor filter		Actions → Tremor filter	Ctrl + F	
Change frequency of tremor filter		Actions → Tremor filter frequency		Select frequency from drop down list
Close pop-up window			ESC	
Change settings in Real Time ECG	Pop-up real time	View → Real Time ECG		Select settings in drop down list
Change settings in Rhythm ECG	Pop-up in rhythm	View → Rhythm ECG		Select settings in drop down list
Change trends in trend display	Pop-up in trends	View → Trends		Select trends in drop down list
Show/hide trend information				Displayed trends depend on <i>Settings</i>








CARDIOLEX

Vretenvägen 6
171 54 Solna
Sweden

Version 4.0

Commands in Report mode

The summary below shows the available functions in the Report mode.

Function	Icon	Menu	Fn	Comment
Change examination data		File → Patient Info	F2	Not possible to change patient data
Print full report		File → Print	F4	
Preview, change printer, report, pages		File → Print...		Temporary change in settings
Preview the report to be printed out		File → Print Preview	Ctrl+W	
Print a page of the currently displayed report		Print page		The <i>Report</i> tab must be selected
Print the full currently displayed report		Print sub report		
Print an interval in Rhythm ECG	Pop-up in rhythm	Print from/to...		Print from time lock to clicking point
Print an interval in ECG	Pop-up in ECG	Print from/to...		Print from time lock to clicking point
Change final report	Pop-up in report			Temporary change of report content
Save examination		File → Save	F5	
Open examination		File → Examinations	F7	
Exit and return to real time mode		File → Real time	F6	
Exit this program		File → Exit		

CARDIOLEX

Vretenvägen 6
171 54 Solna
Sweden

Version 4.0