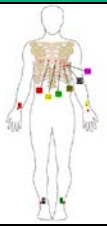

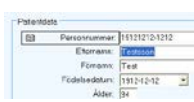
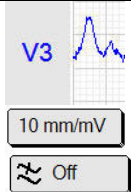







Quick start guide Long ECG Base *EC Sense*

This quick start guide gives a brief explanation on how to use *EC Sense* for acquiring a Long ECG.

Iccon	Fn	Instruction
		1. Prepare the patient and connect all electrodes according to the desired lead electrode placement. Check that the patient cable is connected to the acquisition unit.
	F2	2. Click on the icon for entering patient data.
		3. Enter the patient information in the dialogue box. If the PC is connected to a central storage system the information may be filled in automatically.
		4. Select desired leads, gain and filter settings. Check the quality of the ECG signals.
	F9	5. Click on the red circular icon to start recording an ECG (up to 1 hour can be recorded).
	F9	6. Stop the recording by clicking on the red square icon. Disconnect the patient.
		7. Scroll through the report when it is displayed by clicking on the arrows under the <i>Report</i> tab.
	F4	8. Click on the printer icon to printout the entire recording (or use <i>Print page</i> to print specific pages).
	F5	9. Save the recording on the local PC or in the central storage system.



Read the Operator's manual for complete instructions!

CARDIOLEX

Vretenvägen 6
171 54 Solna
Sweden

Version 4.0