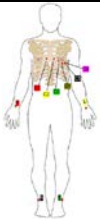

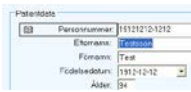
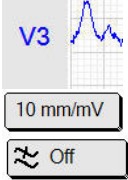
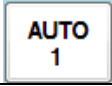
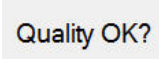




Quick start guide *EC Sense*

This quick start guide gives a brief explanation on how to use *EC Sense* for Resting ECG.

Icon	Key	Instruction
		1. Prepare the patient and connect all electrodes according to desired lead electrode placement. Check that the patient cable is connected to the acquisition unit.
	F2	2. Click on the icon for entering patient data.
		3. Enter the patient information in the dialogue box. If the PC is connected to a central storage system the information will be filled in automatically.
		4. Select desired leads, gain and filter settings. Check the quality of the ECG signals.
	F3	5. Start an automatic sequence, e.g. AUTO 1.
		6. Verify that the quality of the analysed ECG is OK.
	F4	7. Print out the result if desired (when a printer is connected and if not automatically done in the AUTO sequence).
	F5	8. Review and save the analysed ECG on the local PC or in a central storage system. This could be included in the automatic sequence.
		9. Disconnect the patient.



The user should read the entire Operator's manual.

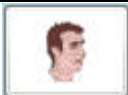

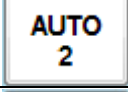








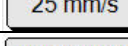
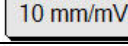


CARDIOLEX

Vretenvägen 6
SE-171 54 Solna
Sweden

Version 4.0

ECG Commands *EC Sense*

The table below shows the commands that are available for Resting ECG.

Function	Icon	Menu item	Fn key	Comment
Enter patient data		<i>File → Patient Info</i>	F2	
Automatic sequence 1		<i>File → Auto 1</i>	F3	
Automatic sequence 2		<i>File → Auto 2</i>	F4	
Print out real time ECG		<i>File → Man Print</i>	F5	
Print out rhythm ECG		<i>File → Rhythm Print</i>	F6	
Stop print out		<i>File → Man Print/Rhythm Print</i>	F5, F6	The button is active if a printout has started
Examinations		<i>File → Examinations</i>	F7	
Settings		<i>File → Settings</i>	Ctrl+ Shift+ Up	Configure the program
Exit this program		<i>File → Exit</i>		Click on X in the top right hand corner
Exercise Test		<i>File → Exercise Test</i>	F8	Only active if Exercise ECG option is installed
Start Long ECG			F9	Only active if Long ECG option is installed
Stop Long ECG			F9	
Leads		<i>View → Leads</i>		Click inside the leads area
Speed		<i>View → Speed</i>		
Gain		<i>View → Gain</i>		
Tremor filter		<i>View → Tremor Filter</i>	Ctrl + F	
mV-pulse		<i>View → Millivolt</i>	F10	
Patients		<i>Tools → Patients</i>		








CARDIOLEX

Vretenvägen 6
SE-171 54 Solna
Sweden

Version 4.0

ECG report commands *EC Sense*

The table below shows the commands that are available in the report mode for Resting ECG.

Function	Icon	Menu item	Fn key
Edit patient information		<i>File → Patient Info</i>	F2
Print out default report		<i>File → Print</i>	F4
Change printout temporarily		<i>File → Print...</i>	
Change report temporarily		Right hand click inside the report for a pop-up menu	
Print the page that is on the screen		<i>→Print page</i>	
Print the report shown		<i>→Print sub report</i>	
Preview the report shown		<i>File → Print Preview</i>	Ctrl + W
Save examination		<i>File → Save</i>	F5
Exit and return to real time mode		<i>File → Real time</i>	F6
Open an examination		<i>File → Examinations</i>	F7
Configure the program		<i>File → Settings</i>	Ctrl+Shift+Up
Exit this program		<i>File → Exit</i>	

CARDIOLEX

Vretenvägen 6
SE-171 54 Solna
Sweden

Version 4.0